

Hiview Plus Microscope

Dual-Usage : Wi-Fi & Wired USB

Adjustable focus and Magnification

5X~200X

ONE WiFi CONNECTION:
WiFi Microscope + Internet

Wi-Fi
1.3MPixel



Leading Innovation>>>

Explore the micro world!
Have fun sharing images!

Professional software with Calibration and Measurement

WiFi Microscope
And Video Camera

The Microscope YOU NEED

Model: MV130WU

Quick Guide: Working at WiFi mode


- Download and install "Hiview Plus" app.
- Slide the Power switch to the "WiFi". Wait about 20 seconds.
- Click WiFi receiving device's wireless connection , select SSID - "WiFi2SCOPE_xxxxxx" network --> connect -->enter password "12345678".
- Open "Hiview Plus" app to see your images.

OS supporting: Windows,Mac,IOS and Andriod

Package includes: WiFi/USB Microscope, USB mini cable, USB AC charger
CD software, Quick Guide, Calibration Ruler,Plastic stand



WiFi Microscope and Video Camera

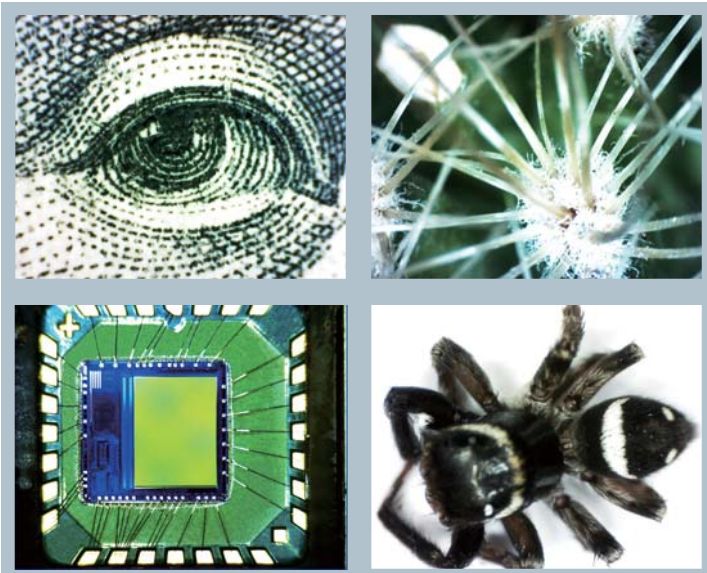

Features:

- Magnification: 5x - 200x @ 17" monitor full screen view
- Optical Format: 1/4" CMOS Sensor(4:3)
- Field of View:
 - Min.: 1.9 mm x 1.5 mm @17 mm working distance(200X)
 - Max.: 82mm x 65 mm @300 mm working distance(5X)
- Image Sensor: 1.3 Mega pixels
- Video Capture Resolution 1280x1024, 640x480
- Image Color: YUV or MJPG
- Lens: High Precision Microscope Lens
- Focus Range: Manual from 8mm to 300mm
- Light Source: Integrated adjustable LED illumination
- Dimension: 142 mm(L) X 36 mm(D) , Weight: 82g

WiFi Specifications:

- IEEE 802.11 b/g/n standards compliant
- Wireless Transmission Distance: up to 5m
- Frequency Range: 2.400 ~ 2.4835GHz
- Lithium Battery: Continuous Working Time - Appx.2 Hours, Full Charging Time - 3 Hours
- Power Source: DC5.0V/1A with mini USB cable

